


VITAMIN D TESTING SIMPLIFIED

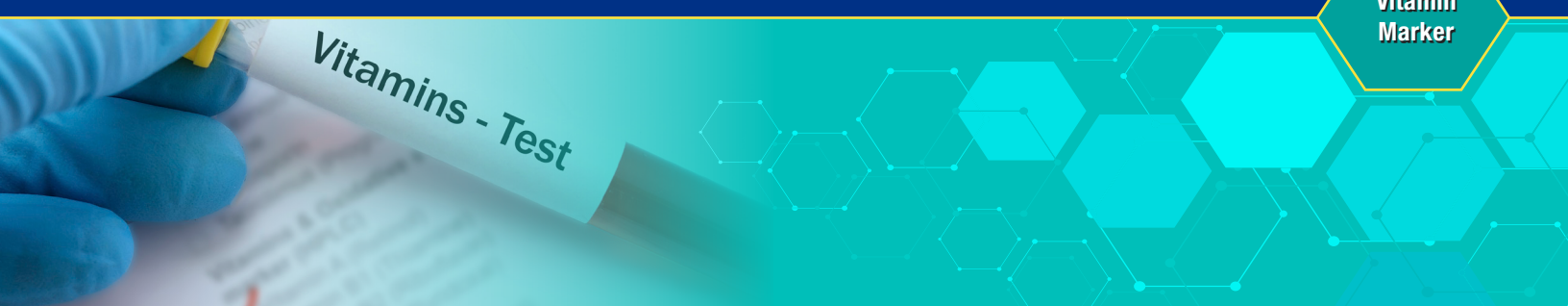

Vitamin
Marker

The *First* and *Only*
TWO REAGENT VITAMIN D ASSAY
FOR CLINICAL CHEMISTRY ANALYZERS

- Recognizes both Vitamin D2 & D3 equally
- Rapid throughput with an assay time <10 min, and reaches more than 500 tests/hr on the Roche Cobas c701
- No sample pre-treatments or pre-dilution steps required
- Low sample volume of 3 μ L
- Good correlation to LC-MS/MS method
- No risk of biotin interference
- Certified by the CDC VDSCP Program
- Meets the performance target set by the DEQAS advisory panel
- Traces to NIST
- Excellent assay precision (%CV <5% at 30 ng/mL)
- Wide assay dynamic range (7.6 – 147.8 ng/mL)



510(k) Cleared; EU:  
2797
Health Canada Registered

25-OH VITAMIN D ASSAY


For Clinical Chemistry Analyzers

Dual Vial
Liquid Stable

ASSAY SPECIFICATIONS

Method	Latex Enhanced Immunoturbidimetric
Sample Type & Volume	- Serum - Plasma Sample Volume 3 µL
Method Correlation	N = 171 y-intercept = -3.03 Slope = 1.062 R² = 0.9575 Samples Ranged From: 8.4 - 146.8 ng/mL
Linearity	7.6 - 147.8 ng/mL
LOB	1.2 ng/mL
LOD	2.9 ng/mL
LOQ	7.6 ng/mL
Calibration Levels	5-Point Calibration

Diazyme EZ Vitamin D Assay Procedure*

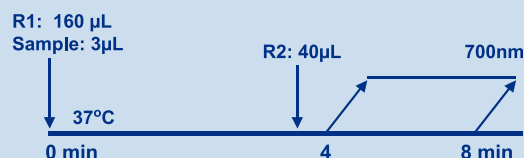

*Analyzer Dependent


For a list of validated parameters please contact Diazyme technical support at 858-455-4768 or email support@diazyme.com

ASSAY PRECISION


25-OH Vitamin D (ng/mL)			Within-Run		Between-Run		Total	
Specimen	n	Mean	SD	%CV	SD	%CV	SD	%CV
Control #1	80	21.7	0.9	3.9%	0.6	2.8%	1.3	6.2%
Control #2	80	42.5	1.0	2.4%	0.8	2.0%	1.7	3.9%
Sample #1	80	22.1	0.8	3.8%	0.8	3.8%	1.2	5.6%
Sample #2	80	42.8	0.9	2.0%	1.0	2.4%	1.3	3.1%
Sample #3	80	99.5	1.8	1.8%	1.5	1.6%	2.7	2.8%

ASSAY METHOD COMPARISON

A comparison of the EZ Vitamin D assay (y) using samples measured with LC-MS/MS (x) gave the following correlation:

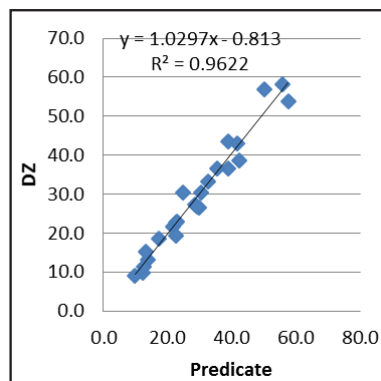


A comparison of the EZ Vitamin D assay (y) using samples measured with a commercially available 25-hydroxyvitamin D immunoassay (x) gave the following correlation:

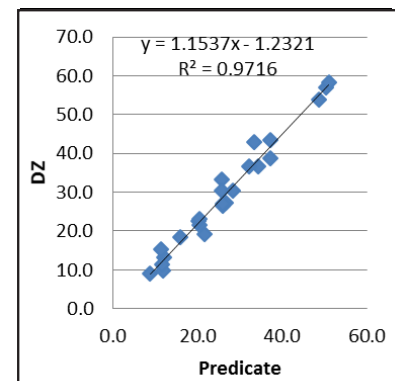

Conclusion: Diazyme EZ Vitamin D assay has a good correlation with a commercially available 25-hydroxyvitamin D immunoassay (R² > 0.97)

VITAMIN D DEFICIENCY & INTOXICATION

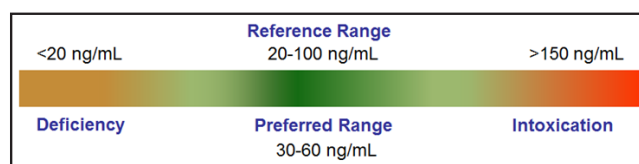

DIAZYME LABORATORIES, INC.

12889 Gregg Court, Poway, CA 92064 USA
PO Box 85608, San Diego, CA 92186 USA
Tel: +1-858-455-4768 +1-888-DIAZYME
www.diazyme.com sales@diazyme.com

DIAZYME EUROPE GMBH

Zum Windkanal 21, 01109 Dresden, Germany
Tel: +49-351-886-3300 Fax: +49-351-886-3366
sales@diazyme.de

SHANGHAI DIAZYME CO., LTD.

Building 8, 3879 Dongchuan Road
Minhang District, Shanghai, 200245
Tel: +86-21-54843802 Fax: +86-21-51320663
www.lanyuanbio.com service@lanyuanbio.com

